

SIGNAL PROCESSING IN MNE: DAY 1

1 – TIMESERIES DATA THE RAW CLASS

07.02.2024

Thomas S. Binns

thomas-samuel.binns@charite.de




EINSTEIN
CENTER
Neurosciences



Brain and
Data Science
QAI LABS




Timeseries data

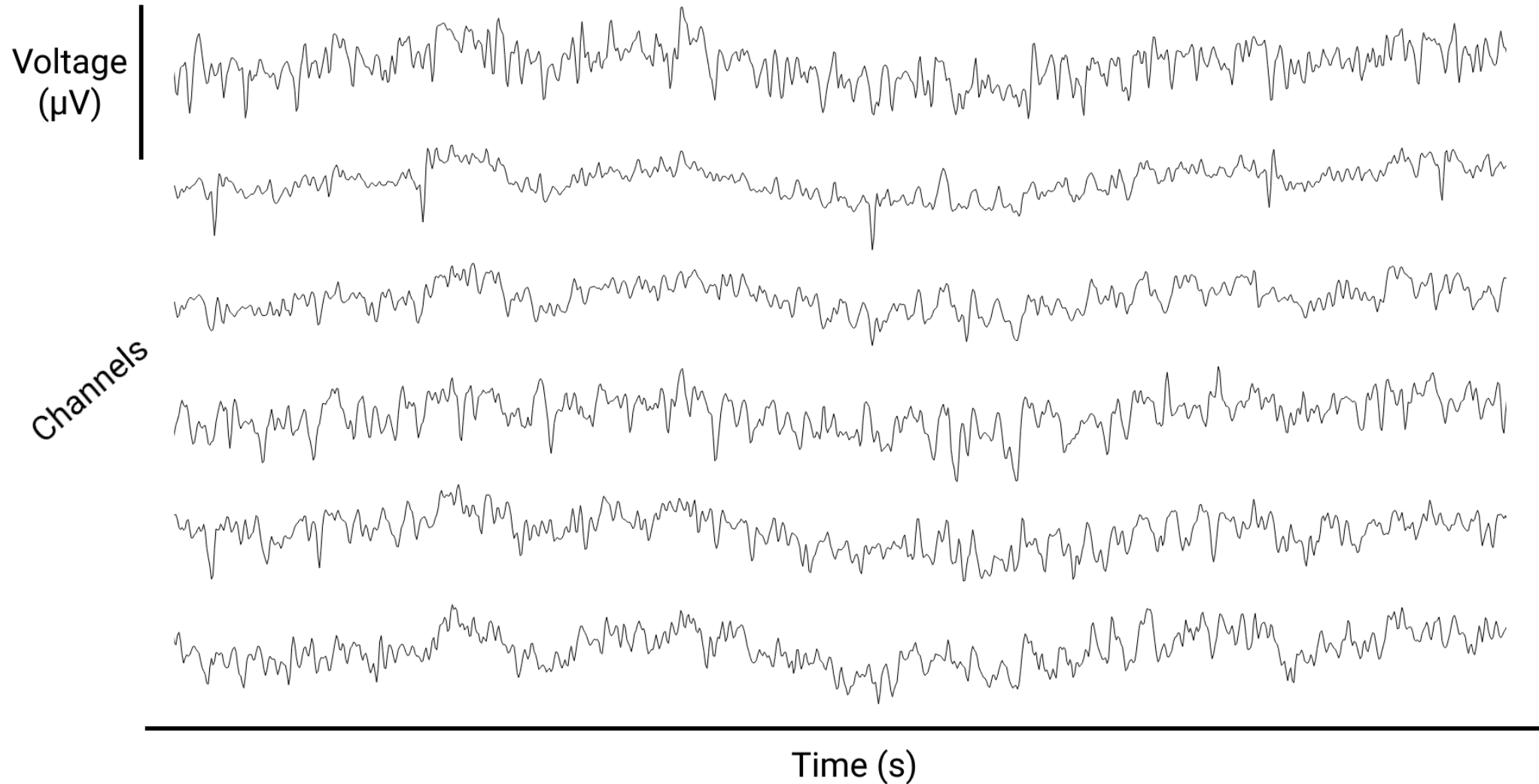

Timeseries data

[Install](#) [Documentation](#) [API Reference](#) [Get help](#) [Development](#)



Section Navigation

- Most-used classes
 - mne.io.Raw**
 - mne.Epochs
 - mne.Evoked
 - mne.Info
- Reading raw data
- File I/O
- Creating data objects from arrays
- Exporting
- Datasets
- Visualization
- Preprocessing
- Events
- Sensor Space Data
- Covariance computation
- MRI Processing

mne.io.Raw

```
class mne.io.Raw(fname, allow_maxshield=False, preload=False,
on_split_missing='raise', verbose=None)
```

[\[source \]](#)

Raw data in FIF format.

Parameters:

fname : *path-like* | *file-like*

The raw filename to load. For files that have automatically been split, the split part will be automatically loaded. Filenames not ending with `raw.fif`, `raw_sss.fif`, `raw_tsss.fif`, `_meg.fif`, `_eeg.fif`, or `_ieeg.fif` (with or without an optional additional `.gz` extension) will generate a warning. If a file-like object is provided, preloading must be used.

Changed in version 0.18: Support for file-like objects.

allow_maxshield : *bool* | *str* (default `False`)

If True, allow loading of data that has been recorded with internal active compensation (MaxShield). Data recorded with MaxShield should generally not be loaded directly, but should first be processed using SSS/tSSS to remove the compensation signals that may also affect brain activity. Can also be "yes" to load without eliciting a warning.

Onto the notebook...

Conclusion

- Timeseries data stored in `Raw` & `RawArray` objects
- Can be loaded from data on disk
`read_raw_xxx()` → `Raw`
- Can be created from data arrays
`array (channels, times)` → `RawArray`
- Various methods for manipulating and visualising data stored in `Raw` objects